

Visualization of NIfTI Images

Formats of Images

There are multiple imaging formats. We will use NIfTI:

- ▶ NIfTI - Neuroimaging Informatics Technology Initiative (<https://nifti.nimh.nih.gov/nifti-1>)
 - ▶ essentially a header and data (binary format)
 - ▶ will have extension .nii (uncompressed) or .nii.gz (compressed)
 - ▶ we will primarily use compressed NIfTI images
 - ▶ we will use 3-D images (4-D and 5-D are possible)
- ▶ ANALYZE 7.5 was a precursor to NIfTI
 - ▶ had a `hdr` file (header) and `img` file (data)

An Imaging into to R

What is a package?

A package is collection of functions, documentation, data, and tutorials (called vignettes).

- ▶ You install a package using the `install.packages` command/function:

```
install.packages("oro.nifti")
```

`install.packages` is a function, "`"oro.nifti"`" is a character string.

Loading Packages

When you install a package, that means it's downloaded on your computer. That **doesn't** mean that you can use the functions from that package just yet.

- ▶ You “load”/import a package into memory using the `library` command

For example, to load the `oro.nifti` package:

```
library(oro.nifti)
```

Now, functions from the `oro.nifti` package can be used.

Some packages we will use

All packages we will discuss are loaded on the RStudio Server:

- ▶ `oro.nifti` - reading/writing NIfTI images
- ▶ `neurobase` - extends `oro.nifti` and provides helpful imaging functions

Let's load them:

```
library(oro.nifti)  
library(neurobase)
```

Reading in NIfTI images: assignment

We will use the `readnii` function (from `neurobase`) to read in a `nifti` object (this is an R object).

Here we read in the `training01_01_mprage.nii.gz` file, and assign it to an object called `t1`

```
t1 = readnii("training01_01_mprage.nii.gz")
```

Now, an object `t1` is in memory/the workspace.

Reading in NIfTI images: assignment

In R, you can assign using the equals = or arrow <- (aka assignment operator).

The above command is equivalent to:

```
t1 <- readnii("training01_01_mprage.nii.gz")
```

There are no differences in these 2 commands, but just personal preference.

By default, if you simply pass the object, it is printed, we can also do `print(t1)`:

```
class(t1)
```

```
[1] "nifti"  
attr(,"package")  
[1] "oro.nifti"
```

```
t1
```

NIfTI-1 format

```
Type : nifti  
Data Type : 4 (INT16)  
Bits per Pixel : 16  
Slice Code : 0 (Unknown)  
Intent Code : 0 (None)  
Qform Code : 1 (Scanner_Anat)  
Sform Code : 1 (Scanner_Anat)  
Dimension : 256 x 256 x 120  
Pixel Dimension : 0.88 0.88 1.17
```

Help

- ▶ To see the documentation for a function, use the ? symbol before the name of the function. This is a shortcut for the help command:
- ▶ For example, to see documentation for readnii:

```
?readnii  
help(topic = "readnii")
```

- ▶ To search for help files, use a double ?? or help.search:

```
??readnii  
help.search(pattern = "readnii")
```

Some Details

- ▶ R is case sensitive (e.g. y and Y are different)
- ▶ Commands separated by new line or by a colon (;
- ▶ Use # to comment

You can also be explicit about which package you are using with the :: operator, where the syntax is package::function:

```
t1 = neurobase::readnii("training01_01_mprage.nii.gz")
```

Multiple Sclerosis

- ▶ Multiple sclerosis is a chronic disease of the central nervous system (brain, spinal chord, optic nerves)
- ▶ Multiple sclerosis lesions in the brain are areas of active inflammation, demylenation, or permanent tissue damage.
 - ▶ lesions are primarily in the white matter
- ▶ MRI is well-suited for assessing lesion burden (volume and patterns) because lesions appear as hyperintensities on FLAIR, T2-w, and PD images and as hypointensities on T1-w images.

`ms.lesion` Package

The `ms.lesion` package was made for this course. It has all the data we will work with and the outputs from the analyses.

The main function we will use is

`get_image_filenames_list_by_subject`. It returns a list of filenames of the images for the 5 training and 3 test subjects from the MS lesion challenge 2016

(<http://iacl.ece.jhu.edu/index.php/MSChallenge>)

- ▶ The full data is available at
<https://smart-stats-tools.org>.

ms.lesion Package

```
library(ms.lesion)
files = get_image_filenames_list_by_subject()
length(files); names(files);
```

```
[1] 8
```

```
[1] "test01"      "test02"      "test03"      "training01"   "tr
[6] "training03"  "training04"  "training05"
```

```
head(files$training01)
```

```
"library/ms.lesion/extdata/raw/training01/training01_01_mpr...
```

```
"library/ms.lesion/extdata/raw/training01/training01_01_mpr...
```

```
"library/ms.lesion/extdata/raw/training01/training01_01_f...
```

Reading in a T1 image

Let's read in the T1 image from a MS lesion data set to visualize:

```
files = files$training01  
t1_fname = files["MPRAGE"]  
t1 = readnii(t1_fname)  
t1
```


NIfTI-1 format

```
Type : nifti  
Data Type : 4 (INT16)  
Bits per Pixel : 16  
Slice Code : 0 (Unknown)  
Intent Code : 0 (None)  
Qform Code : 1 (Scanner_Anat)  
Sform Code : 1 (Scanner_Anat)  
Dimension : 256 x 256 x 120  
Pixel Dimension : 0.83 x 0.83 x 1.17  
Voxel Units : mm  
Time Units : sec
```

Orthographic view

The `oro.nifti::orthographic` function displays nifti objects in 3 different planes.


```
oro.nifti::orthographic(t1)
```

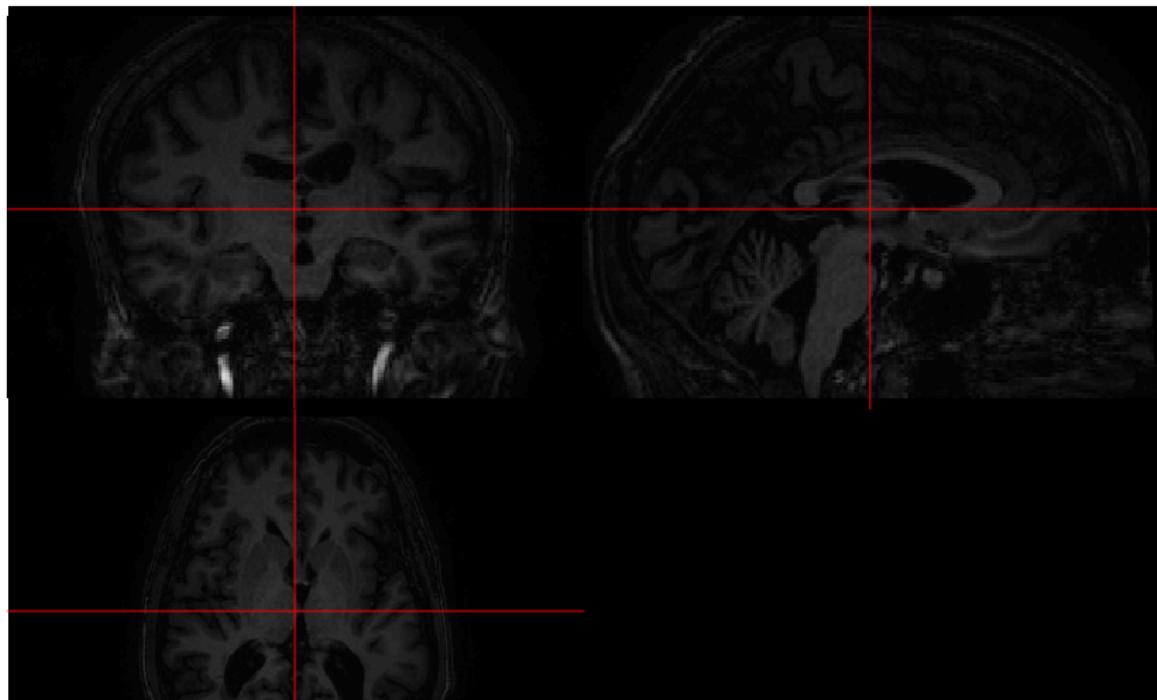

Orthographic view with additions

The `neurobase::ortho2` function expands upon this with some different defaults.


```
neurobase::ortho2(t1)
```

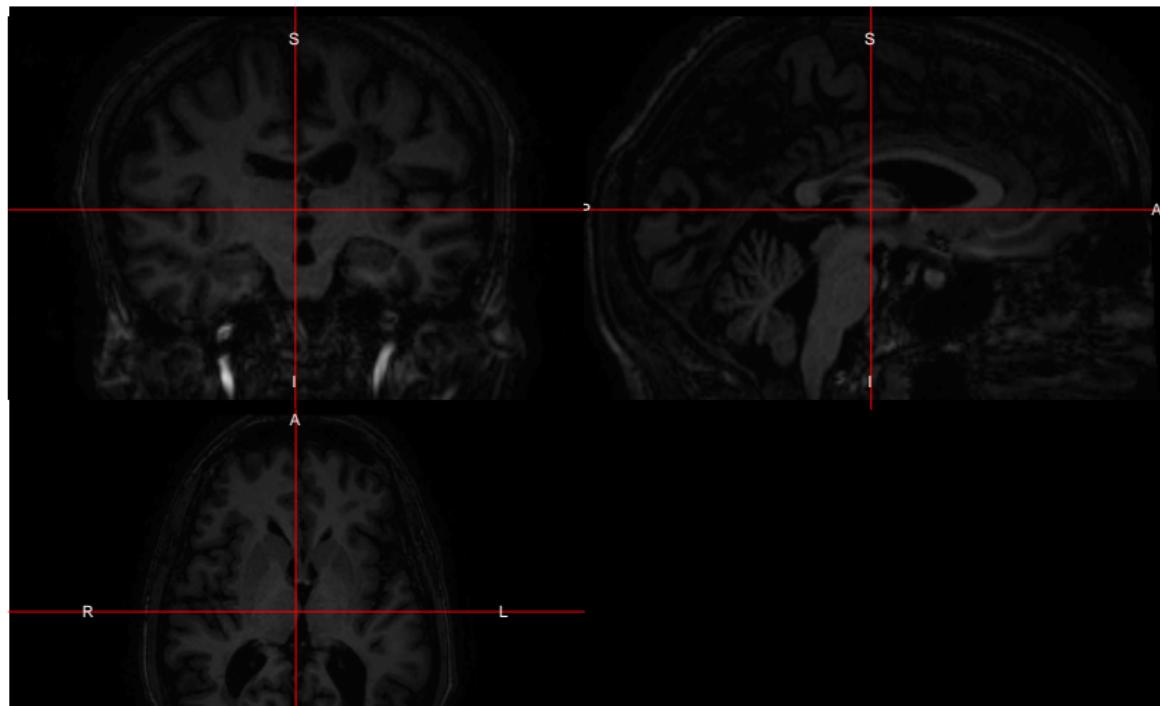

Orthographic view with additions

ortho2 annotates the orientation of the image. These many not be correct. You can turn these off with the add.orient argument:


```
neurobase::ortho2(t1, add.orient = FALSE)
```

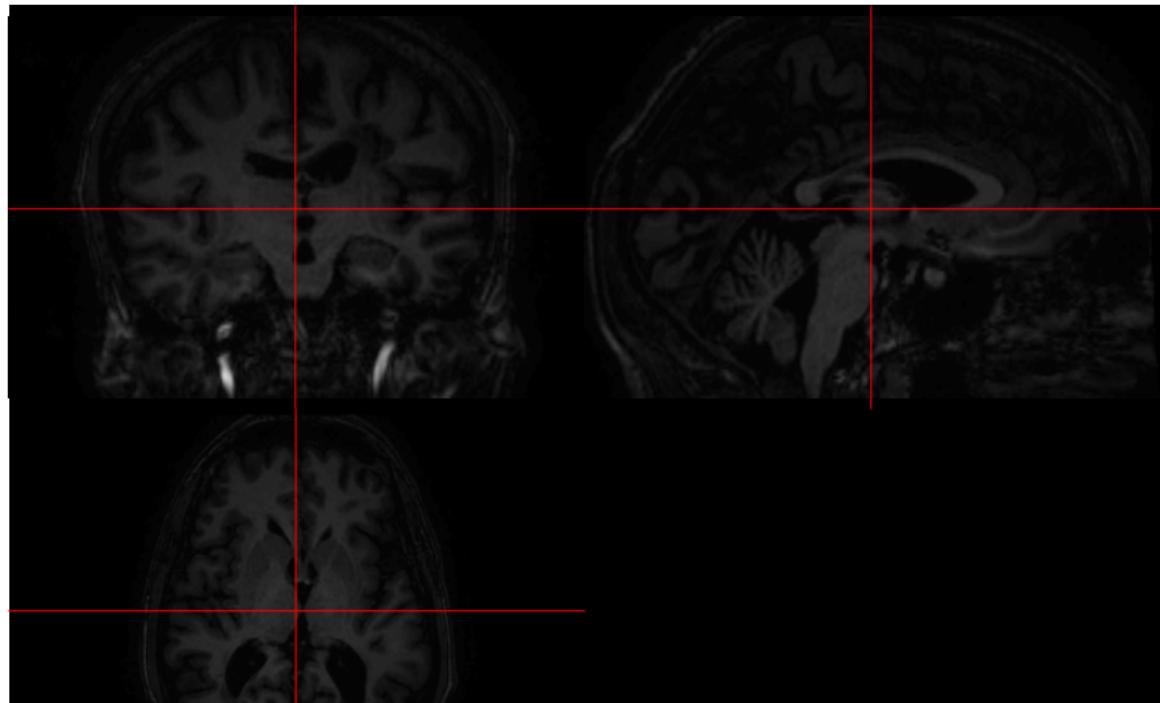

Bright values


High intensities can “dampen” the viewing of an image.

`ortho2(t1)`

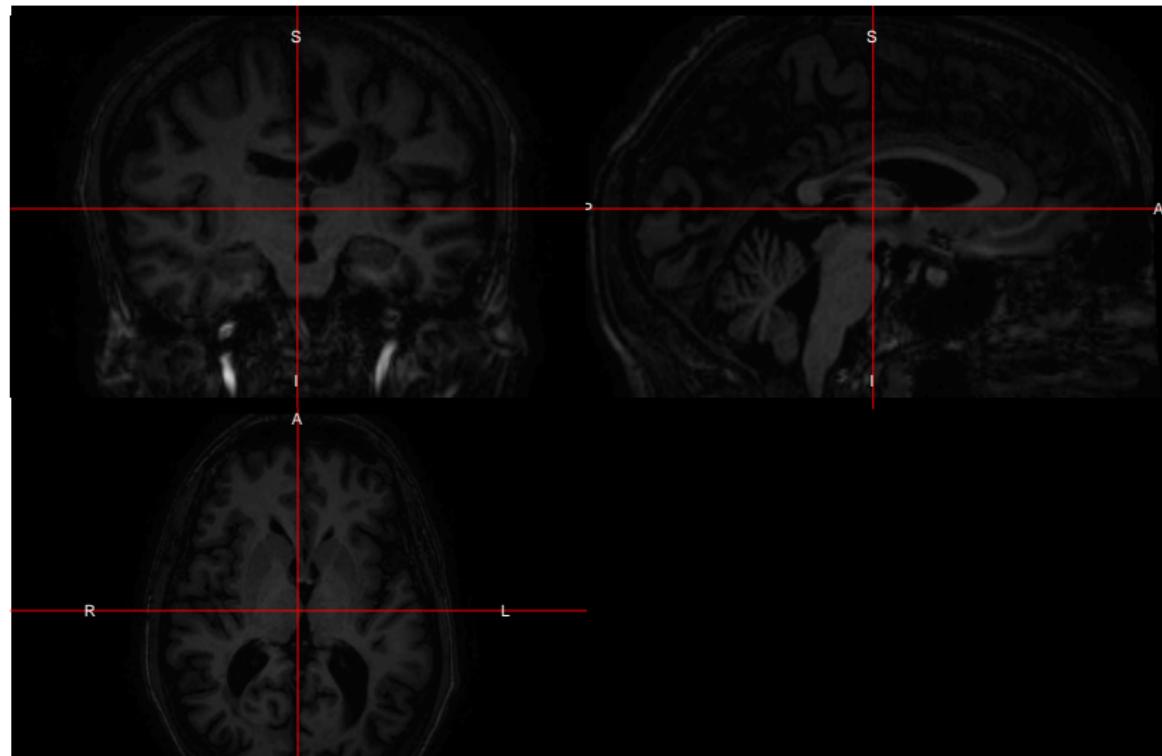

We see a faint outline of the image, but this single large value affects how we view the image. The function `robust_window` calculates quantiles of an image, by default the 0 (min) and 99.9th quantile, and sets values outside of this range to that quantile.


```
ortho2(robust_window(t1))
```

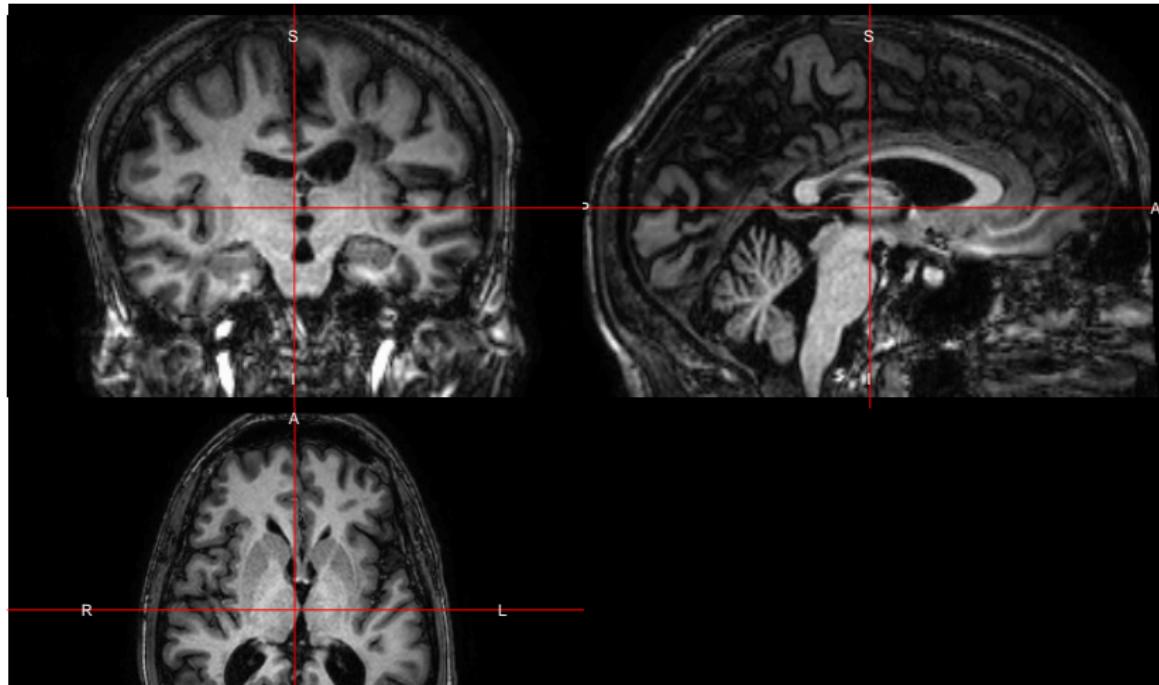

Changing the Windowing

The `zlim` option can also map intensities that can be plotted:

```
ortho2(t1, zlim = quantile(t1, probs = c(0, 0.999)))
```

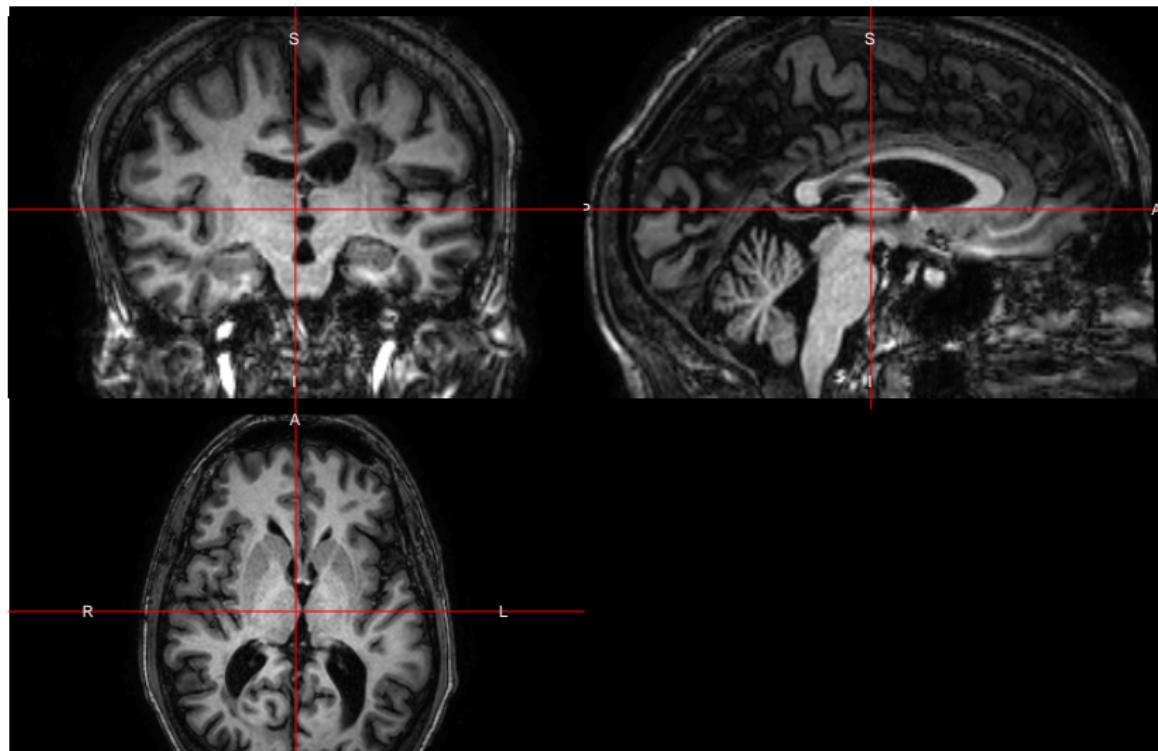

Differences between orthographic and ortho2

The above code does not fully illustrate the differences between `orthographic` and `ortho2`. One marked difference is when you would like to “overlay” an image on top of another in an orthographic view.

We will be plotting the windowed image for visualization:

```
rt1 = robust_window(t1)
```


Here we will highlight voxels greater than the 90th quantile of the image:

```
qmask = t1 > quantile(t1, 0.9)
```

Differences between orthographic and ortho2

We see that the white matter is represented here, but we would like to see areas of the brain that are not over this quantile to be shown as the image.


```
orthographic(rt1, y = qmask)
```

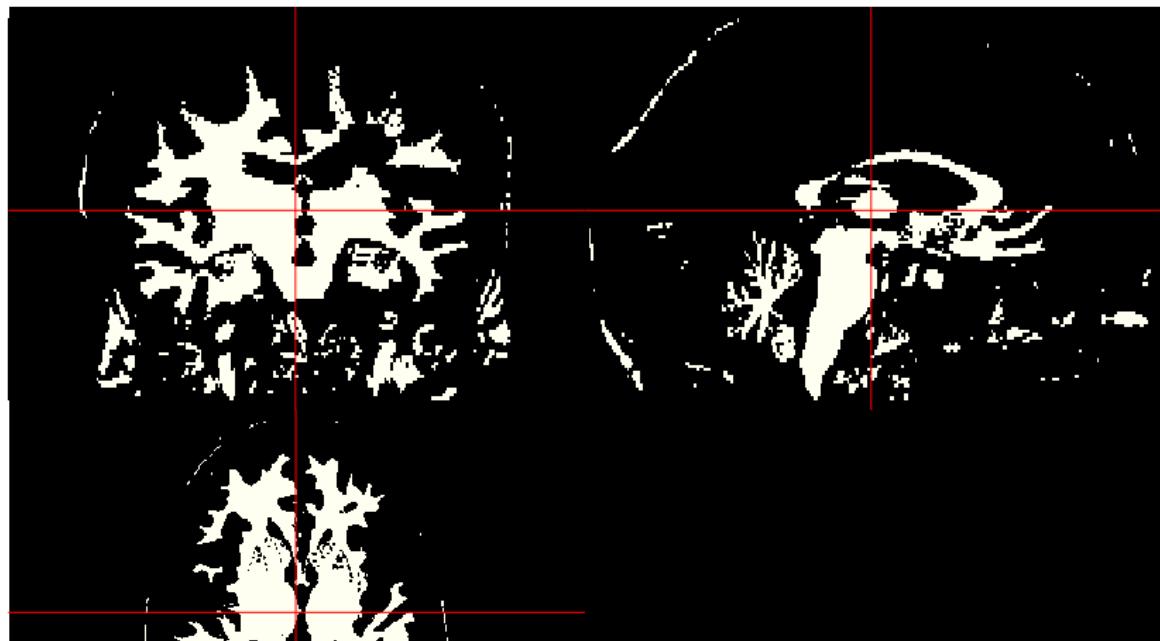

Differences between orthographic and ortho2

Let us contrast this with:

```
ortho2(rt1, y = qmask)
```

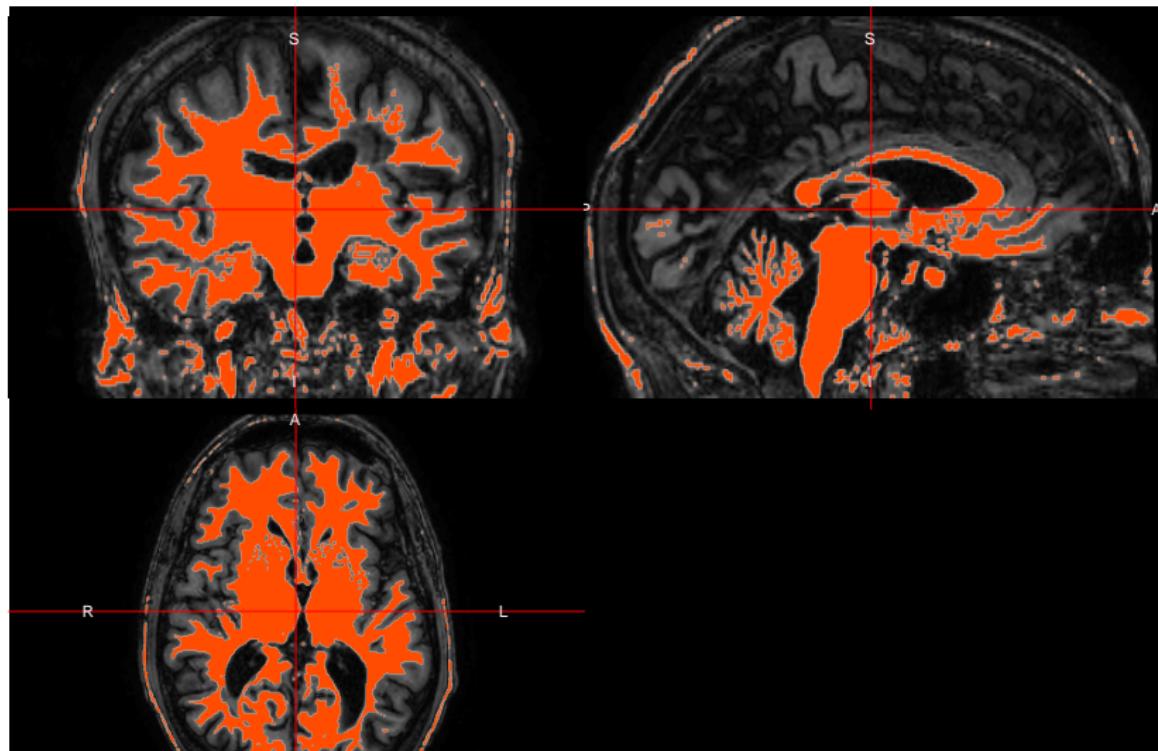

Differences between orthographic and ortho2

In R, NA is the placeholder for missing data.

The ortho2 (and orthographic) function is based on the graphics::image function in R, as well as many other functions we will discuss below. When graphics::image sees an NA value, it does not plot anything there.

We see the image where the mask is 0 shows the original image. This is due to the NA.y argument in ortho2. The NA.y argument in ortho2 makes it so any values in the y argument (in this case the mask) that are equal to zero are turned to NA.

Double orthographic view


Sometimes you would like to represent 2 images side by side, of the same dimensions and orientation of course. The `double_ortho` function allows you to do this. Let's get a coarse mask of the image using the `oMask` function:

```
mask = extrantsr::oMask(t1)
```

Double viewing images

We can view the original T1 alongside the estimated mask:

```
double_ortho(t1, mask)
```

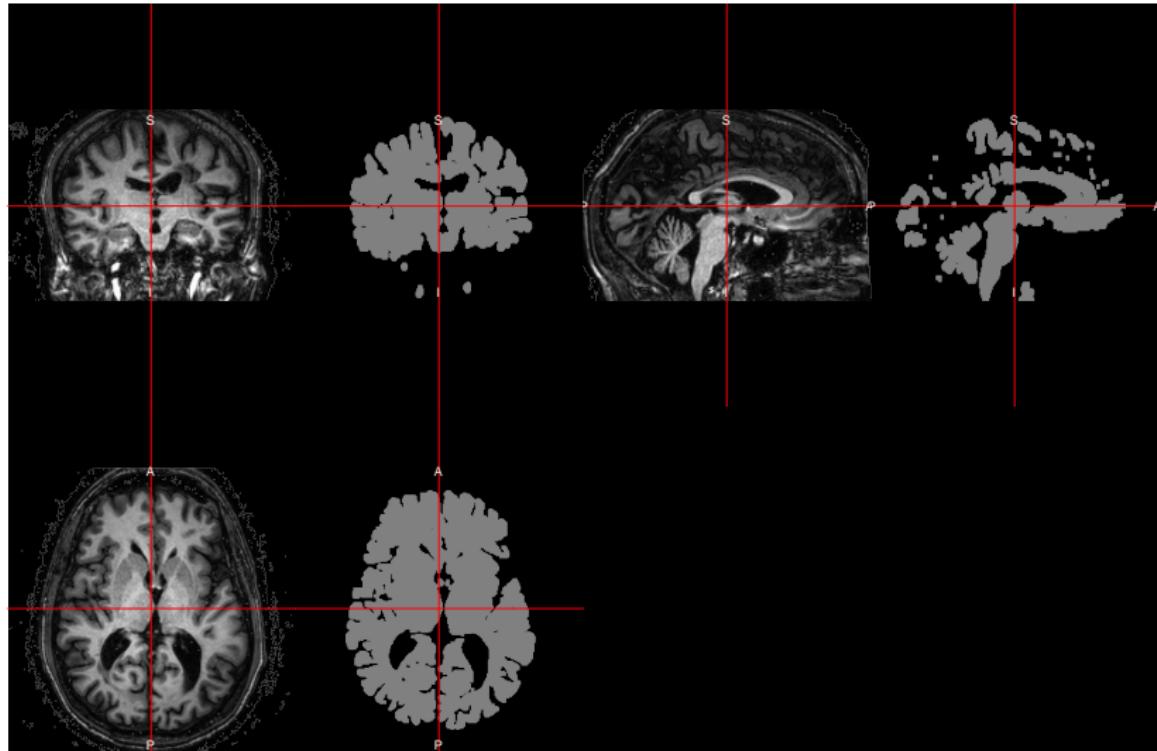

Single slice view

We may want to view a single slice of an image. The `oro.nifti::image` function can be used here. Note, `graphics::image` exists and `oro.nifti::image` both exist.


The `oro.nifti::image` allows you to just write `image(nifti_object)`, which performs operations and calls functions using `graphics::image`.

This allows the user to use a “generic” version of `image`, which `oro.nifti` adapted specifically for `nifti` objects. You can see the help for this function in `?image.nifti`.

Plotter

Let's **try** to plot an image of the 80th slice of the T1 image:


```
image(t1, z = 80)
```

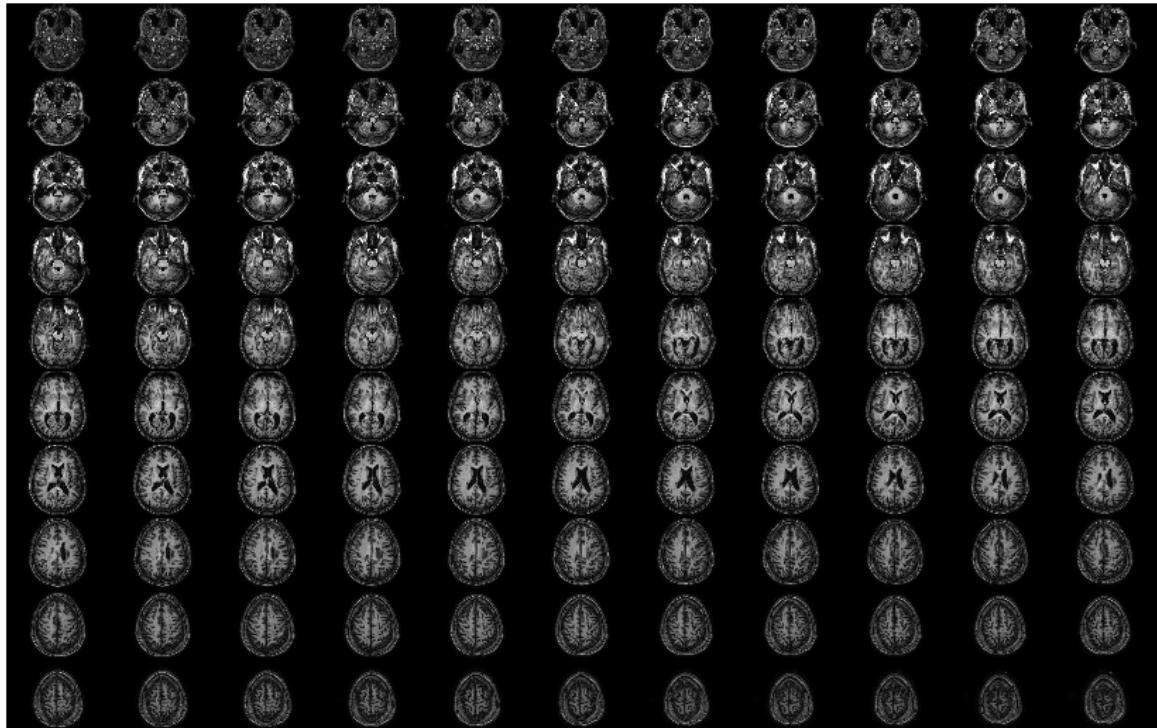

Bad Plot

What happened? Well, the default argument `plot.type` in `image.nifti` is set for "multiple", so that even if you specify a slice, it will plot **all** slices. Here, if we pass `plot.type = "single"`, we get the single slice we want.


```
image(t1, z = 80, plot.type = "single")
```

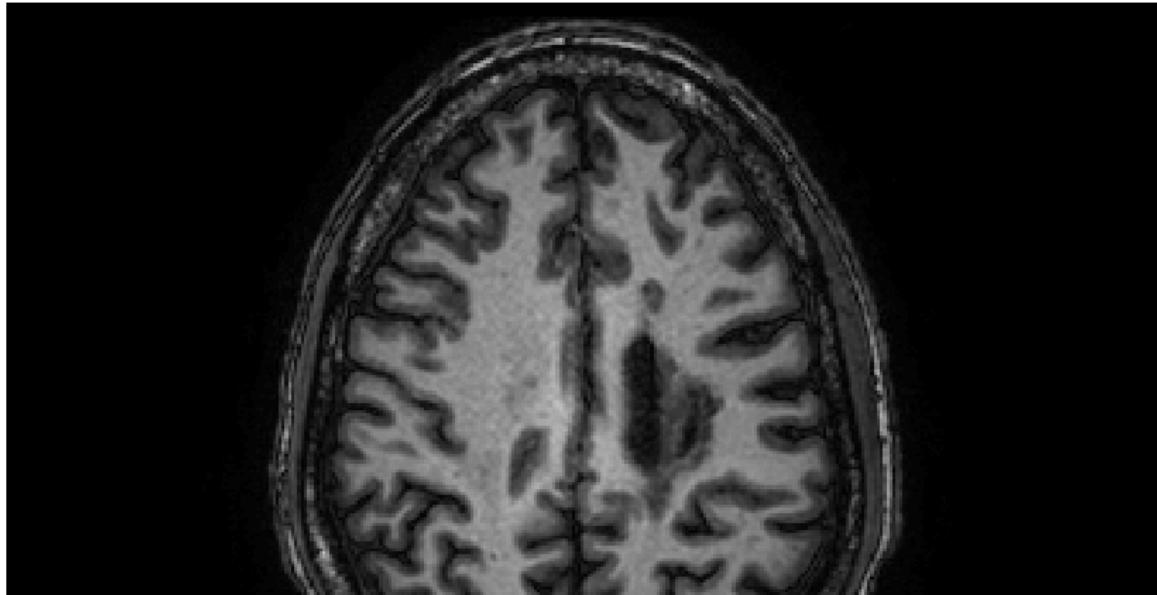

Multiple Slices but not All

If we put multiple slices with `plot.type = "single"`, then we will get a view of these 2 slices:


```
image(t1, z = c(60, 80), plot.type = "single")
```

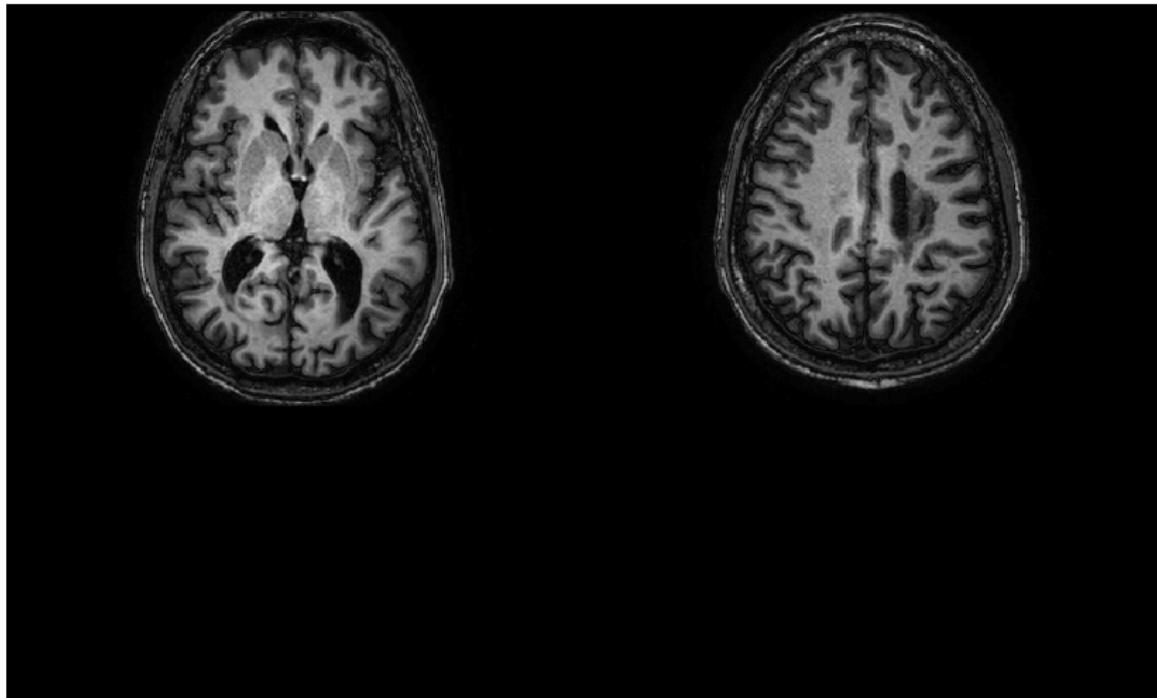

Multiple Slices but not All

In oro.nifti version 0.7.4, we added slice, which is the same as image, but plot.type = "single" as default:

```
oro.nifti::slice(t1, z = c(60, 80))
```

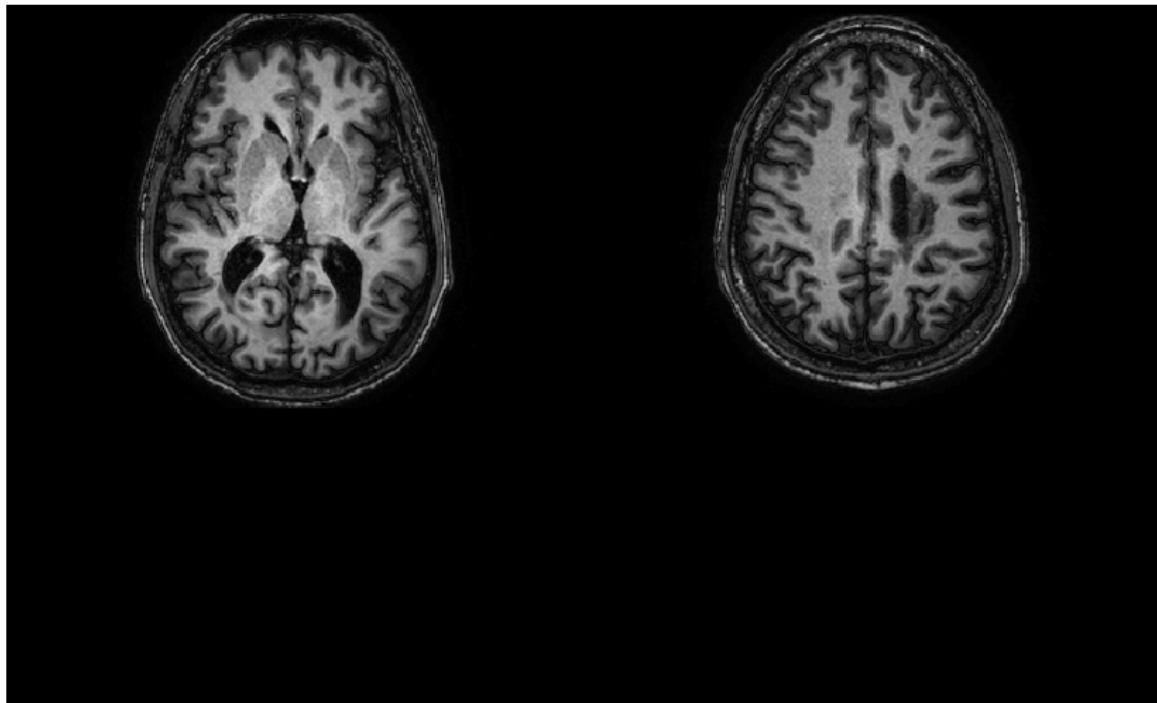

Overlaying slices

We can also overlay one slice of an image upon another using the `oro.nifti::overlay` function. Here we must specify `plot.type` again for only one slice.


```
overlay(t1, y = qmask, z = 80, plot.type = "single")
```



Overlay2

There is no overlay2 because we added NA.y in oro.nifti version 0.7.4:

```
overlay(t1, y = qmask, z = 80, plot.type = "single", NA.y =
```

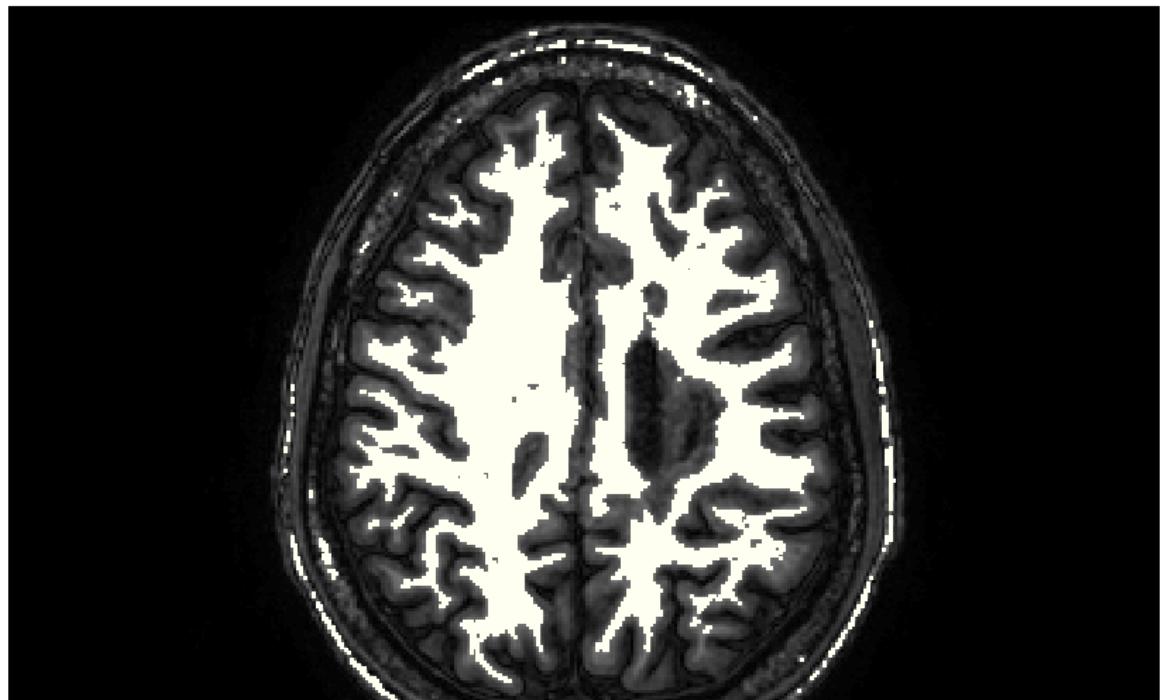

Dropping empty dimensions

In some instances, there are extraneous slices to an image. For example, in the Eve template image we read in, it is just the brain. Areas of the skull and extracranial tissue are removed, but the slices remain so that the brain image and the original image are in the same space with the same dimensions. For plotting or further analyses, we can drop these empty dimensions using the `neurobase::dropEmptyImageDimensions` function or `drop_empty_dim` shorthand function.

By default, if one `nifti` is passed to the function and `keep_ind = FALSE`, then the return is a `nifti` object.

```
reduced = dropEmptyImageDimensions(t1)
dim(t1)
```

```
[1] 256 256 120
```

```
dim(reduced)
```

```
[1] 256 256 120
```

Plotting the Reduced Image

ortho2(reduced)

